

JH-351R4 Digital Programmable USB Rechargeable Hearing

Aid / Hearing Amplifier

Specification of JH-351R BTE RIC FM Rechargeable Hearing Aid Hearing Amplifier			
Max Sound Output	125±3dB	Noise EIN	≤20dB
Sound Gain	40±5dB	Battery size	3.7V 30mAh Built-in lithium batter
Total Harmonic Wave Distortion	500Hz ≤ 3.1%	Working current	≤2.5mA
	800Hz ≤ 1.7%		
	1600Hz ≤ 0.2%		
Frequency Range	200-5000Hz	Full charging Time	2 hours
Dimension		Duration of use	More than 20 hours of use
Certification	FDA	Cycle time:	No decline within 500 times charging

Features of JH-351R BTE RIC FM Rechargeable Hearing Aid Hearing Amplifier

1. 2 hours continuous use with a full charge
2. Super mini in the ear, invisible in the ear.
3. Wide frequency range, high quality sound.
4. Open fit hearing aid, no blockage effect, improving wear comfort.
5. Supplied with USB charging dock
6. Different sized ear tips supplied to ensure a good fit

Precautions of JH-351R BTE RIC FM [Rechargeable Hearing Aid Hearing Amplifier](#)

1. Charging time :It will take around 2 hours for a finished battery to be fully recharged. It's strongly suggest to charge for 12 hours for the first time charging.
2. During the charging period, the indication light in the adapter base is showing Red light
3. After charging completed, the indication light will turns to Green light.
4. Important to off the machine before you put on it to your ear.
5. To choose the right size of ear plugs will be more comfortable to the ear,and help a better sounds gain, and avoid any extra noise.
6. The switch should be turned to the off position when the apparatus is not in use. That will ensure longer battery life

JINGHAO

JINGHAO • TO BE BETTER

Package Included

1 [Hearing Aid](#)


3 Ear tips

1 USB Dock

1 Strong box

1 Manual book

2 Cleaning tools

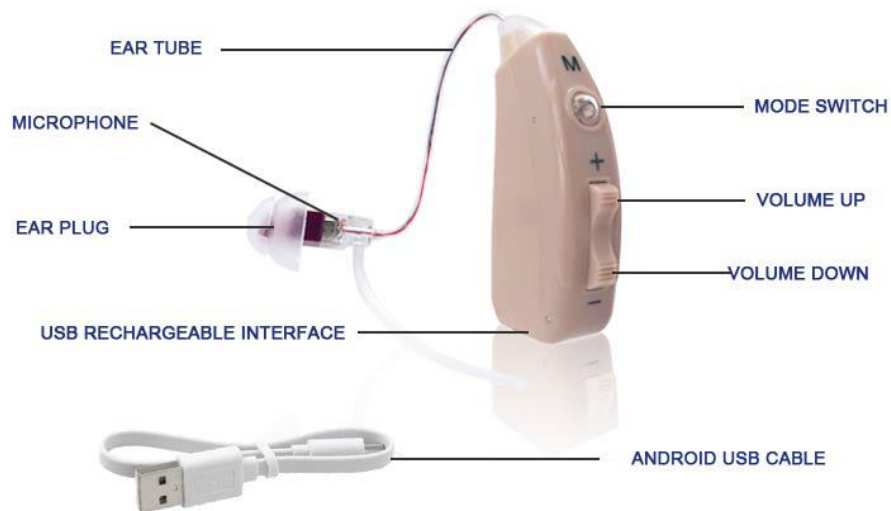

JH-351 RIC

- ✓ Receiver In The Canal
- ✓ High Sound Quality
- ✓ USB Rechargeable
- ✓ Touch-tone Volume Control